
Luxeo 6Z





Let There Be Light!

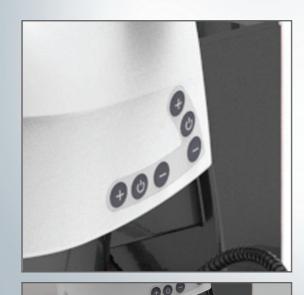

Traditional approaches to stereo illumination can be cumbersome and get in the way of routine operation. Alleviating this pain point inspired a user friendly, versatile approach to delivering brilliant illumination features in the Luxeo 6Z, all in one sleek package.

The Luxeo 6Z allows users to choose an ultra-bright front mounted gooseneck arrangement that efficiently hides all wiring in the arm bracket, while keeping the goosenecks away from focus controls. This provides the user with hassle free control of indirect illumination angles. Ring light only options are also available.

A built-in 23-LED ring light with controllable sectors means users can light up a range of objects with bright, homogeneous, and oblique light. An intuitive light control cluster ensures easy intensity control for both gooseneck and ring light options, situated near the microscope's smooth focus controls for ultimate convenience. As is standard with every Labomed microscope, superb optics are immediately evident when looking through the Luxeo 6Z. Featuring proprietary MaxLite[™] coatings, this microscope offers unparalleled light homogeneity and image accuracy across all magnifications.

A flex-arm system is available to maximize the degree of articulation and control in workbench applications, with two articulating arms and a focus mount. This system features gooseneck and ring light illumination to ensure that your specimen is properly lit, regardless of its topology.

A host of camera accessories are available to make the Luxeo 6Z a documentation machine. Choose from the many Labomed iVu camera options or opt for third party video adapters.







Arm Systems That Will Delight


1 A range of boom stands is available in three different sizes and features to accommodate a rich mix of user needs and budgets. Designed to accommodate small to large samples and parts for inspection, the boom stands are further complemented with focus mounts for more precise optical inspection. All incident illumination features are available. **2** For routine industrial use, a bright field/dark field illuminated base is available, which can be quickly removed and modified into a slim base stand. Users on shop floors can opt for the slim base stand with either a traditional arm or pole stand. Also available is a type sliding mirror BF/DF base for most critical illumination of transparent organisms.

3 The bright field/dark field transmitted illumination base quickly becomes a slim base by detaching the electronic cartridge, so users benefit from two stand types in one. Binocular observation comes standard, with an optional trinocular viewing head style in case digital documentation is required.




Arm and Stand Concepts for an Array of Applications

4 For the most routine of applications, a pole stand with slim base is available for a larger working area to accommodate a range of sample sizes and types. Gooseneck and ring light incident illumination are standard.

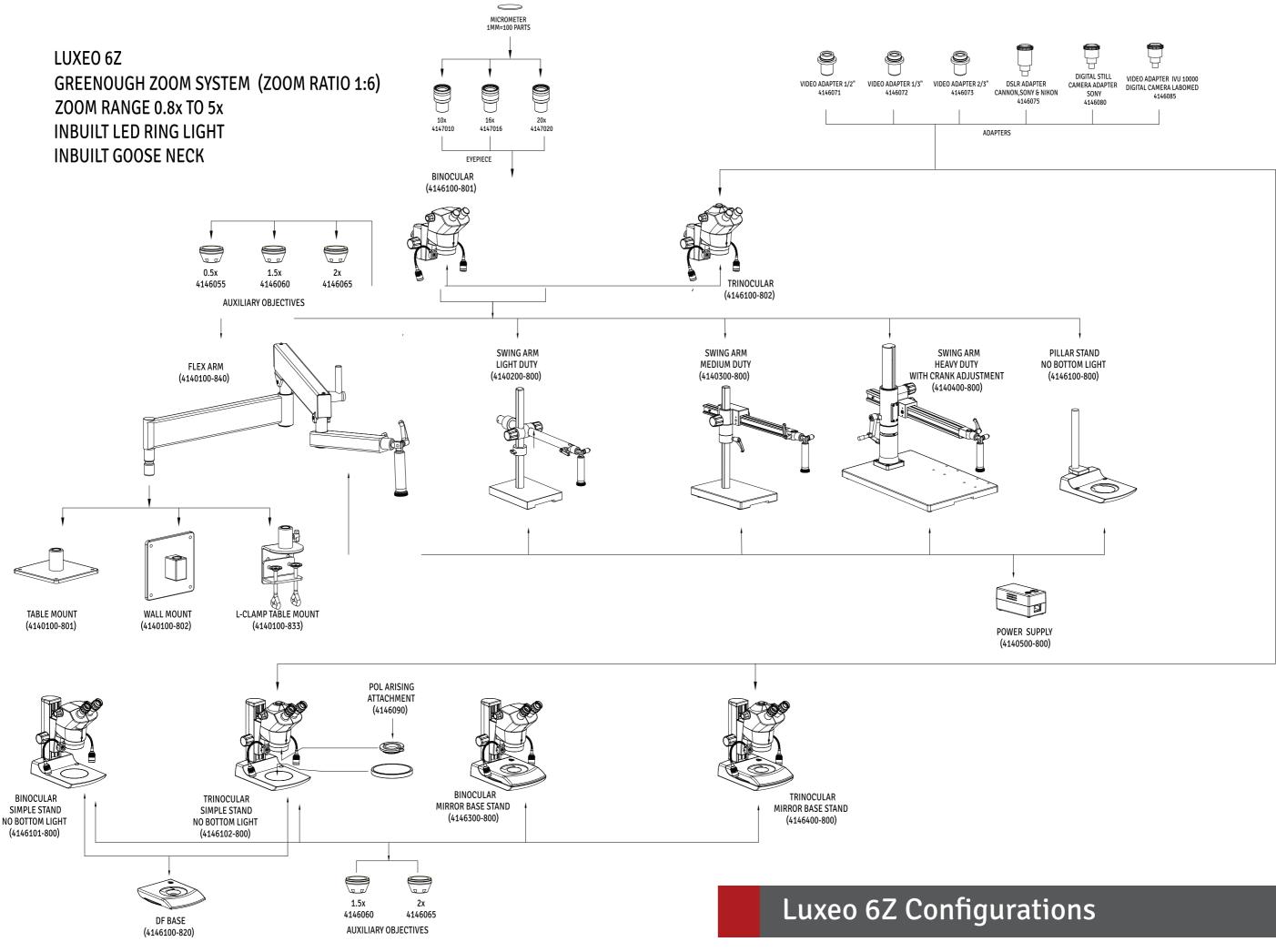

Luxeo 6Z

Optical System Gre	enough
Magnification Zoo	m 6:1
Eyepieces WF	10x/22mm with foldable eye guards, ±5 mm diopter adjustment
Opt	ional: 16X and 20X eyepieces
IPD 50-	75 mm
Viewing Angle Bin	ocular, 45° and Trinocular, 45°
Working Distance 110) mm
Magnification Range 8x t	to 50x
Object Field Diameter 4.4	to 27.5 mm
	dent Light - LED Ring Light with Arc and Circular illumination options
wit	h intensity control:
- Fr	ont Angle Arc illumination with cluster of 11 LEDs
- Re	ear Angle Arc illumination with cluster of 12 LEDs
- Ci	rcular Illumination with cluster of 23 LEDs
- Go	oose neck for extra illumination
Transmitted Light Two	options - detachable BF/DF (Bright field/Dark field) base;
pre	configured mirror base for transparent objects
Auxiliary Objectives Opt	ional: 0.5x, 1.5x, 2.0x
Reticules For	measurement and counting
	ional: CCD adapters 0.35x, 0.50x, 0.65x; DSLR adapters for Canon, Sony and on; Digital camera adapters for Sony and Labomed iVu.
NIK	on, Digital camera adapters for Sony and Laborned IVd.

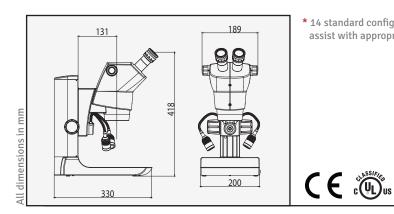

* 14 standard configurations available. Contact a Labomed representative to assist with appropriate configuration selection.

Labo America, Inc. 920 Auburn Ct. Fremont, CA 94538 U.S.A. Tel: (510) 445-1257 Fax: (510) 445-1317 E-mail: sales@laboamerica.com **www.laboamerica.com**



